



4141 MacArthur Blvd. • Newport Beach, CA 92660
 800-854-7256 • Fax 800-411-9722 • glidewelldental.com

Dr. Name _____ Phone # _____

Patient Name _____ Acct. # _____
First Last

Address _____
City/State/ZIP

Email _____ **Deliver by 5 p.m. on** [See Reverse for Working Times](#)

Enclosed with Case: Impressions Models Bite Photos Other:

- Diagnostic Wax-Up**
- Crown
- Veneer No-Prep Veneer
- Additive
- Clear Stent Putty Matrix*
- Open bite _____ mm

SHADE INSTRUCTIONS

Indicate Shade Here

- Transition Crowns & Bridges**
- Transition C&B are milled from an acetal resin block (No reinforcement)

Fabricate:

- Lab Prep Solid Model
- Cement-on Implant #(s) _____
- Screw-retained Implant #(s) _____
- Open bite _____ mm

Smile Transitions

- Upper Lower
- Buy 2, same arch only
- Occlusal holes No occlusal holes
- Open bite _____ mm
- Add length to gingival on tooth
 #(s) _____ by _____ mm

TOOTH NUMBER(S) _____ to _____

Abutment* #(s) _____

Pontic #(s) _____

Veneer #(s) _____

Splinted* #(s) _____ Individual units #(s) _____

Gingival Contour #(s) _____ by _____ mm

Perio treatment: Prepare tooth below gingival on tooth #(s) _____ by _____ mm

Pontic site healing: Prepare ovate socket on tooth #(s) _____ by _____ mm

Amount of Prep Reduction:

- 1 mm 1.5 mm* 2 mm
- Open Trench
- #(s) _____

SELECT INCISAL EDGE TYPE



- Rounded Squared Pointed
- LVI Design Style (DWAX Only)


SELECT PONTIC DESIGN TYPE



- *
- OSDESIGN Image Approval



To expedite your case, please fill out all necessary information.

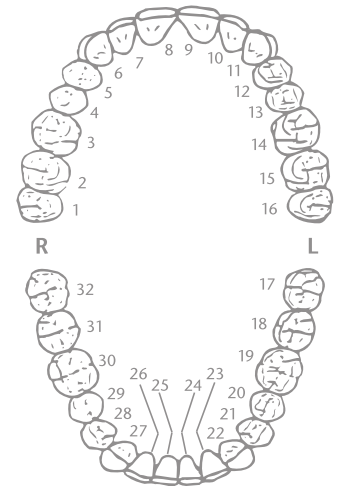

Signature _____

License _____ Date _____

Submission of this Rx constitutes agreement with limited warranty terms and conditions. See reverse for details.

TERMS AND WARRANTY INFORMATION



All Restorations Made in the USA

We honor VISA, MASTERCARD, AMEX and DISCOVER.

TERMS: Cost of collection of any account will be paid by the customer. All accounts are payable within 30 days of statement date. **Accounts not paid within the stated terms will be subject to COD status and a late charge of 2 percent of the unpaid balance.** Prices subject to change without notice. Rx must be enclosed with original case submission.

NO-FAULT REMAKE POLICY: Glidewell is pleased to process all remakes or adjustments at no additional charge if requested within the warranty period and **accompanied by the return of the original appliance.**

LIMITED WARRANTY/LIMITATION OF LIABILITY. For warranty terms and conditions and limitation of liability, visit glidewell.com/policies-and-warranties.


- **Diagnostic Wax-Up (DWAX)**, 1 month warranty
- **Transition Crowns & Bridges**, 2 year warranty
- **Smile Transitions**, 1-year warranty



Please include upper and lower impressions or models and a bite registration. For best results, please do a 'two-step' impression. You may also submit a digital scan. For a detailed process on how to take excellent impressions, scan the QR code to the left.


PREPARATION GUIDELINES

ANTERIOR

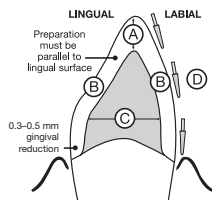


- A. 1 mm incisal reduction
- B. 0.8 mm middle third reduction
- C. Labial and lingual walls must be convergent
- D. Preparation should be cut in three planes

POSTERIOR

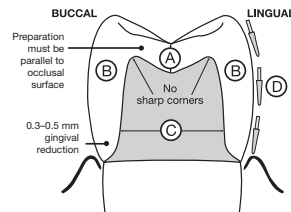


- A. 1 mm occlusal reduction
- B. 0.8 mm middle third reduction
- C. Buccal and lingual walls must be convergent
- D. Preparation should be cut in three planes

LABIAL OR 360° MARGIN

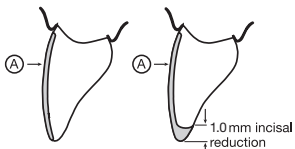

- A. 0.3 to 1.0 mm labial reduction

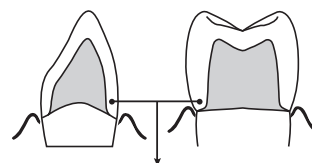

- 0.3-0.5 mm gingival reduction using rounded shoulder margin design


IN-LAB WORKING TIMES

Please allow full working time for each project selected. Working times are **NOT** guaranteed and do **NOT** include weekends or holidays.

Diagnostic Wax-Up (DWAX)	5 working days
Transition Crowns & Bridges	6 working days
over implant	6 working days
Smile Transitions	7 working days


Rush service available. Call for details.

TRANSITION CROWNS & BRIDGES RELINE INSTRUCTIONS

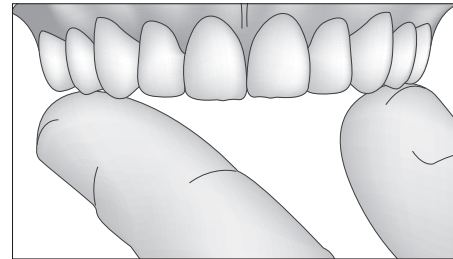


Step 1: If the Transition Crowns & Bridges (TCAB) has already been tried in the mouth, make sure the intaglio is completely dry and free of saliva and debris.

Step 2: Sandblast the intaglio of the TCAB with aluminum oxide. If a benchtop lab sandblaster unit is not available in the practice, a chairside unit can be purchased from dental supply companies such as Henry Schein or Pearson. If sandblasting is not an option, the intaglio of the TCAB can be roughened with a diamond bur. The sandblasting or roughening of the intaglio of this acetal appliance is a necessary step that should not be skipped. Finish with a blast of air or air/water to remove the remaining aluminum oxide particles from the appliance.

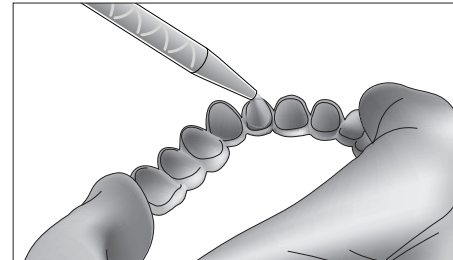


Step 5: Fill the TCAB with bis-acryl material (such as Luxatemp or EXACTA Temp Xtra) and seat the appliance over the preparation(s). Once seated into place, use a microbrush or explorer to remove the excess.

Step 6: Before the bis-acryl material has hardened, light-cure the entire TCAB to adhere the bonding agent/acetyl resin restoration/bis-acryl material sandwich. 5 seconds per unit, per surface is advised.

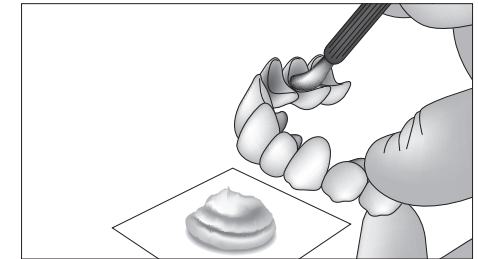

Step 3: Apply a layer of Scotchbond Universal Adhesive (or any other bonding agent containing MDP) to the intaglio of the TCAB. Scrub (and add more bond if necessary) until the intaglio has a slick and shiny look to it (but not dripping wet) and does not evaporate or absorb the bonding agent. Do NOT cure.

Step 4: Block out any undercuts on the actual tooth preparation(s), then apply a layer of Vaseline to the actual tooth preparation(s) in order to avoid adhering the TCAB to the tooth/teeth.

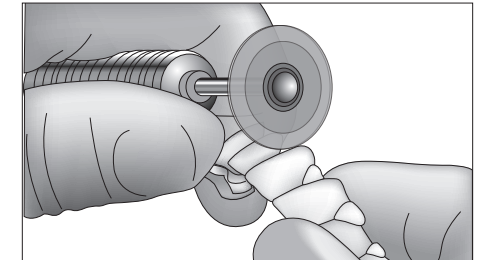

Step 7: Upon removal of the TCAB from the mouth after the bis-acryl material has set, light-cure each unit of the TCAB for an additional 10 seconds through the intaglio surface.

Step 8: Trim and adjust margins and embrasures of the TCAB, check and adjust occlusion, and seat with your provisional cement of choice.